

血清挥发

on MEGA-ALC 1 Column

Column: **MEGA-ALC 1 - 0.20mm, 1.10 μ m, 10m**

Catalog Code: C-ALC1-020-110-10


条件:

进样口: 分流 250 $^{\circ}$ C, 100mL/min分流比, 0.2 μ L进样量

检测器: FID 250 $^{\circ}$ C

柱温箱程序: 35-70 $^{\circ}$ C, 10 $^{\circ}$ C/min

载气: 氦气, 60kPa

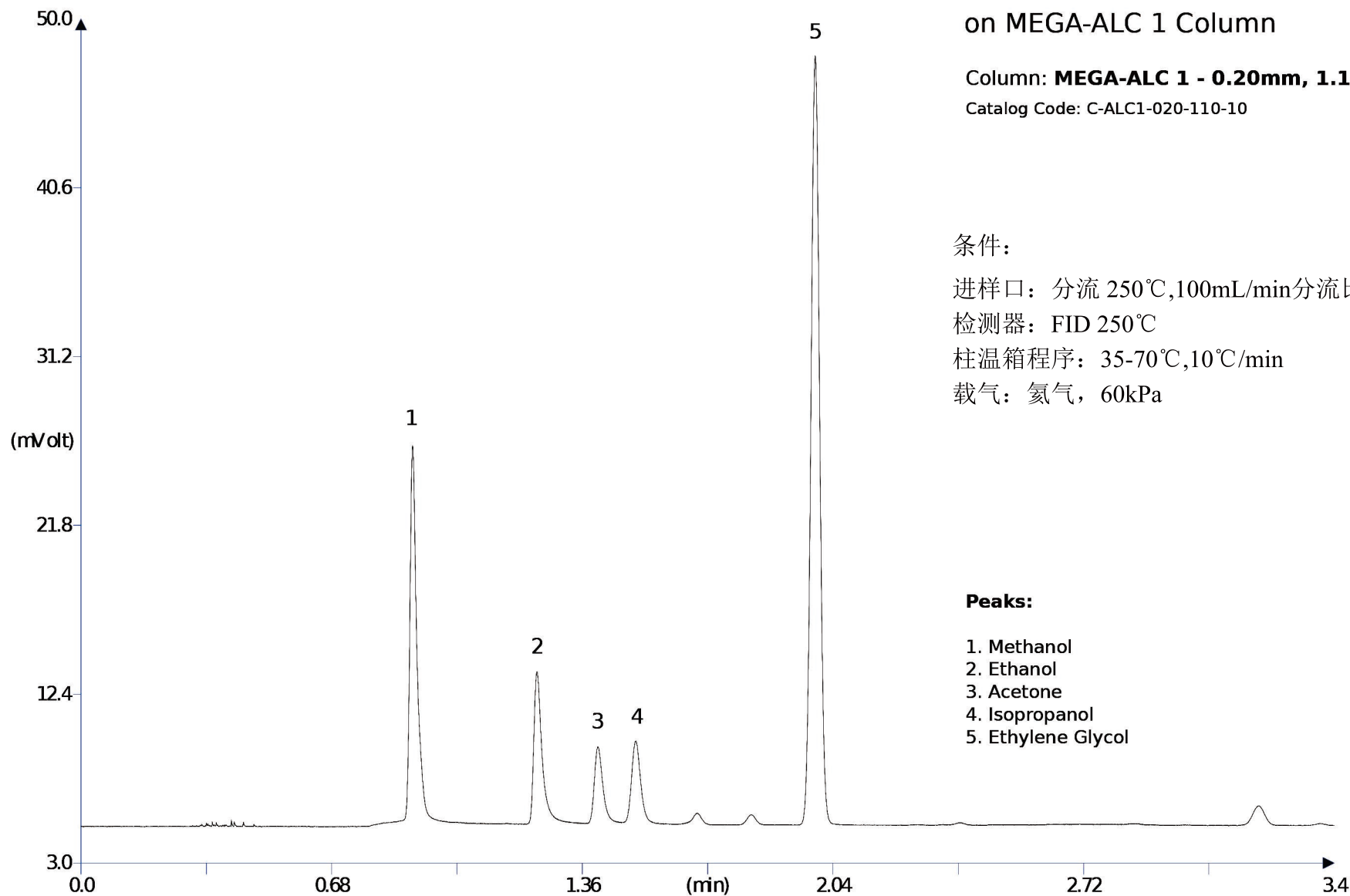

Peaks:

1. Methanol
2. Ethanol
3. Acetone
4. Isopropanol
5. Ethylene Glycol