


(R/S) 1,2-Propanediol - chiral separation

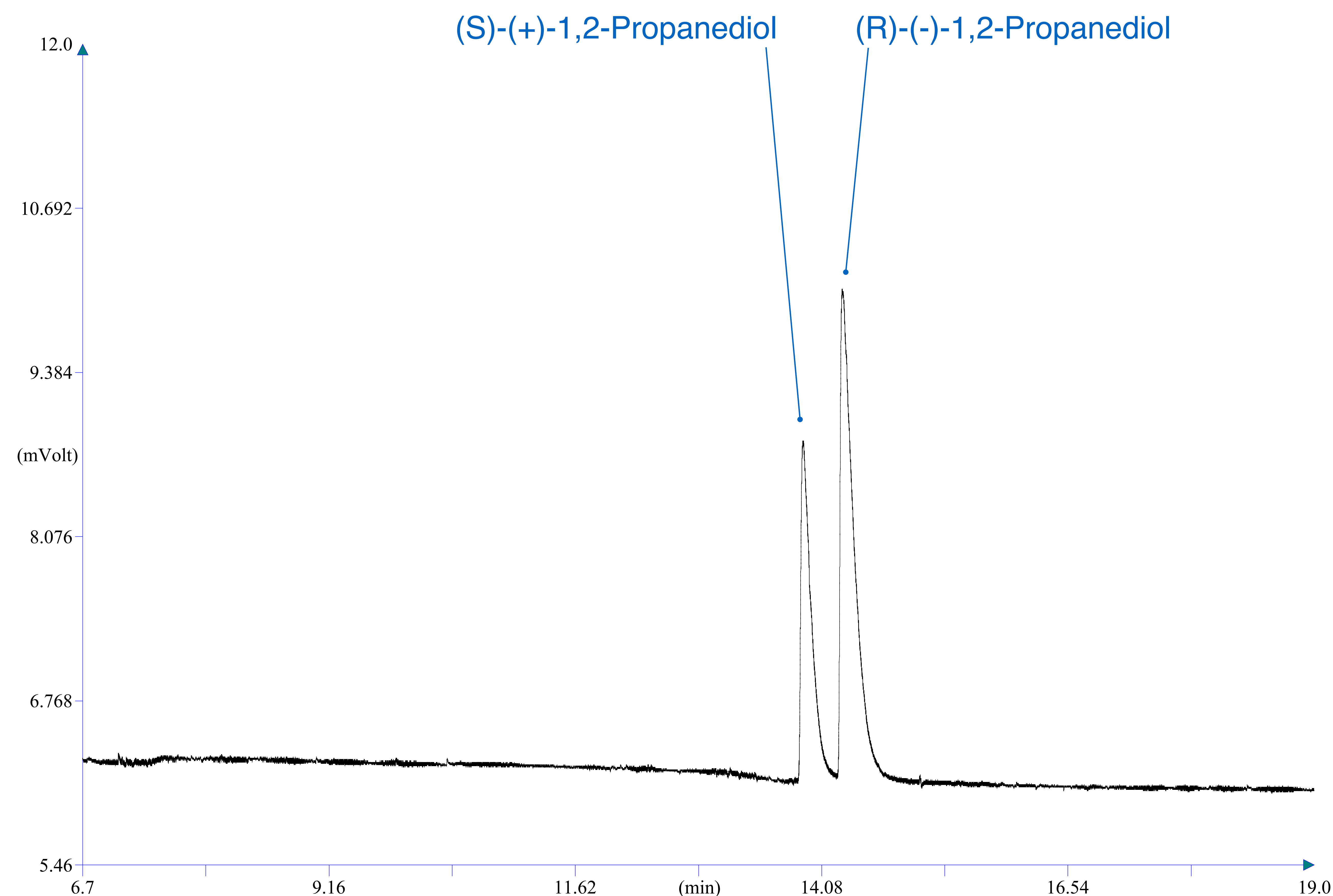

Column

MEGA-DEX DMT - 0.25mm - 0.25micron - 30m
(Catalog Code: DEX-DMT-B-025-025-30)

Conditions

Oven: 50°C, 1.5°C/min.
Carrier: 80 kPa, Hydrogen.
Injector: Split, 250°C.
Detector: FID, 250°C.

Contact us for more details: info@mega.mi.it



MEGA[®]
gas chromatography solutions

2000

MEGA[®]

1980

MEGA[®]

www.mega.mi.it

MEGA S.r.l. - Via Plinio, 29 - 20025 Legnano (MI) - Italy
Tel./Fax: +39 0331.54.79.24 - info@mega.mi.it