


2-苯丙醇异构体手性分离

dex xeb
chiral columns

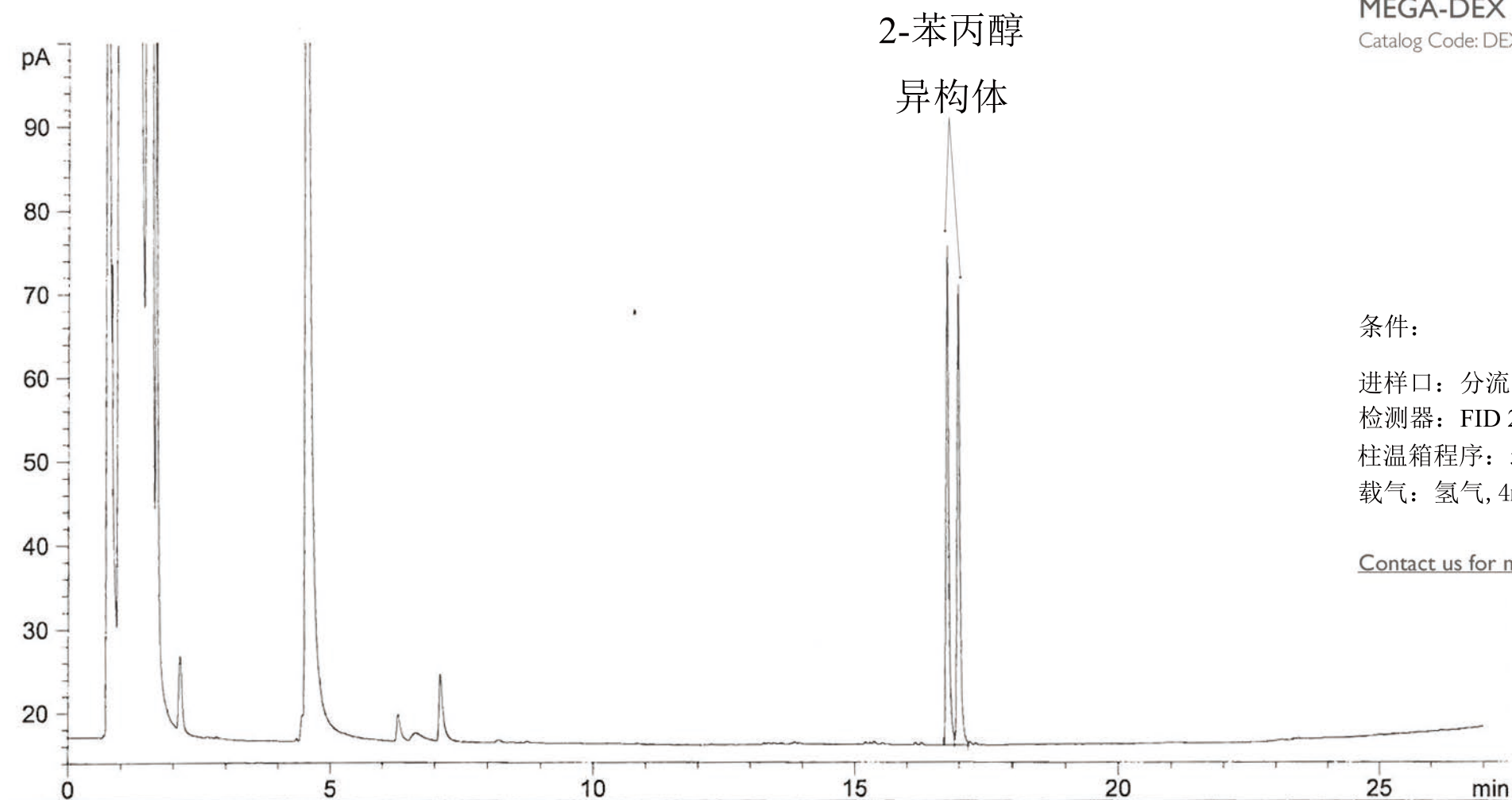

Column:

MEGA-DEX B-SE - 0.25mm, 0.25 μ m, 30m

Catalog Code: DEX-BSE-025-025-30

条件:

进样口: 分流 190 $^{\circ}$ C

检测器: FID 210 $^{\circ}$ C

柱温箱程序: 50 $^{\circ}$ C, 5 $^{\circ}$ C/min

载气: 氢气, 4mL/min

[Contact us for more info](#)

Acknowledgments: author of this work is **Dr. A. Adler**, University of Köln, Organic Chemistry Dept., Greinstrasse 4, Köln, Germany. Many thanks also to and **Dr. R. Forche**, "Dr. R. Forche Chromatographie" (www.gc-hplc.de).

MEGA[®]

www.mega.mi.it