

α -松油醇的手性分离

GC-MS
columns

dex xeb
chiral columns



Column:

MEGA-DEX DET Beta - 0.25mm, 0.25 μ m, 25m

Catalog Code: DEX-DET-B-025-025-25


条件:

进样口: 分流250 $^{\circ}$ C, 50mL/min分流比, 1 μ L的进样量

样品稀释: 1mL溶液中稀释1mg样品

柱温箱程序: 50 $^{\circ}$ C, 2 $^{\circ}$ C/min

载气: 氦气, 1mL/min

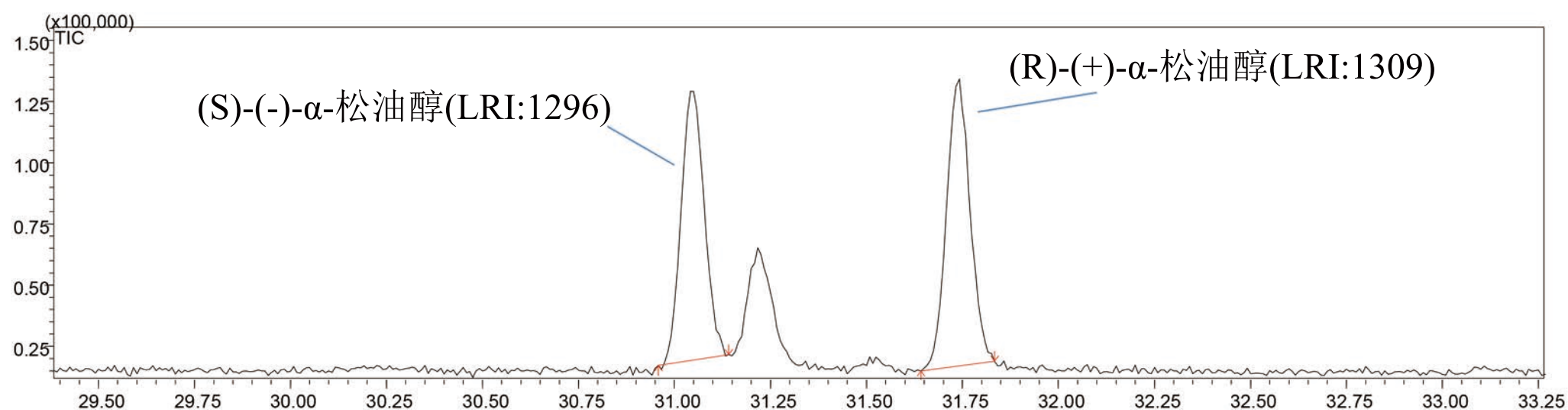

All our MEGA-DEX chiral GC columns work perfectly with MS detectors

Acknowledgments: many thanks to Barbara Sgorbini, Cecilia Cagliero, Chiara Cordero and Carlo Bicchi, Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via P. Giuria, 9, Torino, Italy.

MEGA[®]

www.mega.mi.it