


Column: MEGA-DEX DAC Beta - 0.25mm, 0.25 μ m, 25m
Catalog Code: DEX-DAC-B-025-025-25

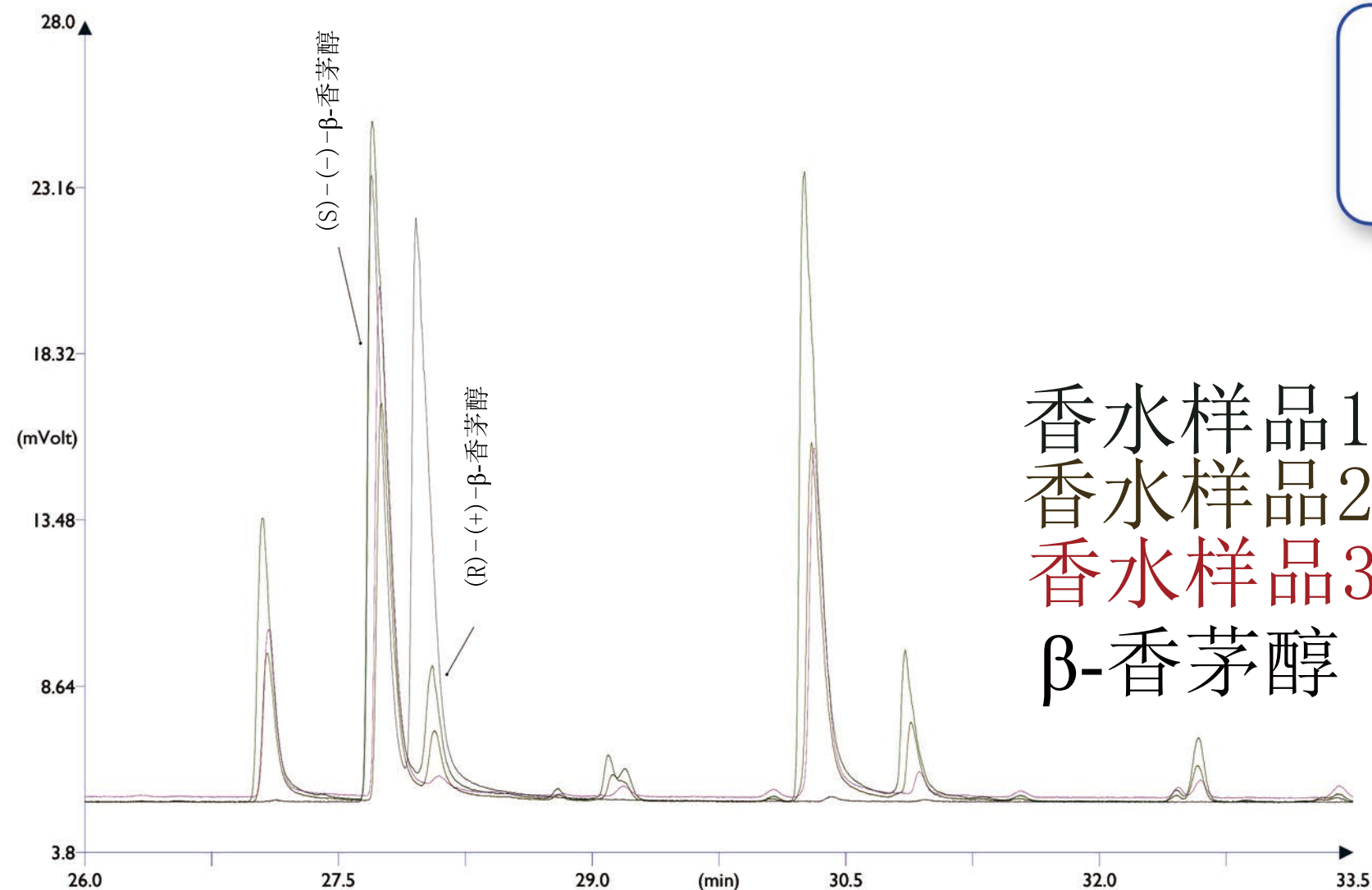

Determination of
(S)-(R)-Beta-Citronellol ratios
in different real perfume samples

香水样品1
香水样品2
香水样品3
 β -香茅醇

条件:

柱温箱程序: 50°C-225°C, 2°C/min

载气: 氢气, 80KPa, 恒压模式

进样口: 分流250°C, 分流比90mL/min, 0.5 μ L的进样量

检测器: FID 250°C