

δ-酯类

Chiral Separation on MEGA-DEX Columns

Column: **MEGA-DEX DAC Beta - 0.25mm, 0.25μm, 25m**

Catalog Code: DEX-DAC-B-025-025-25


条件:

进样口: 分流250℃, 0.5 μL进样量, 1:50分流比

制样: 标准混合, 1mg的样品溶解在1mL的正己烷中

峰:

1. (+/-)δ-辛内酯
2. (+/-)δ-诺娜内酯
3. (+/-)δ-癸内酯
4. (+/-)δ-十一烷酸内酯
5. (+/-)δ-十二内酯

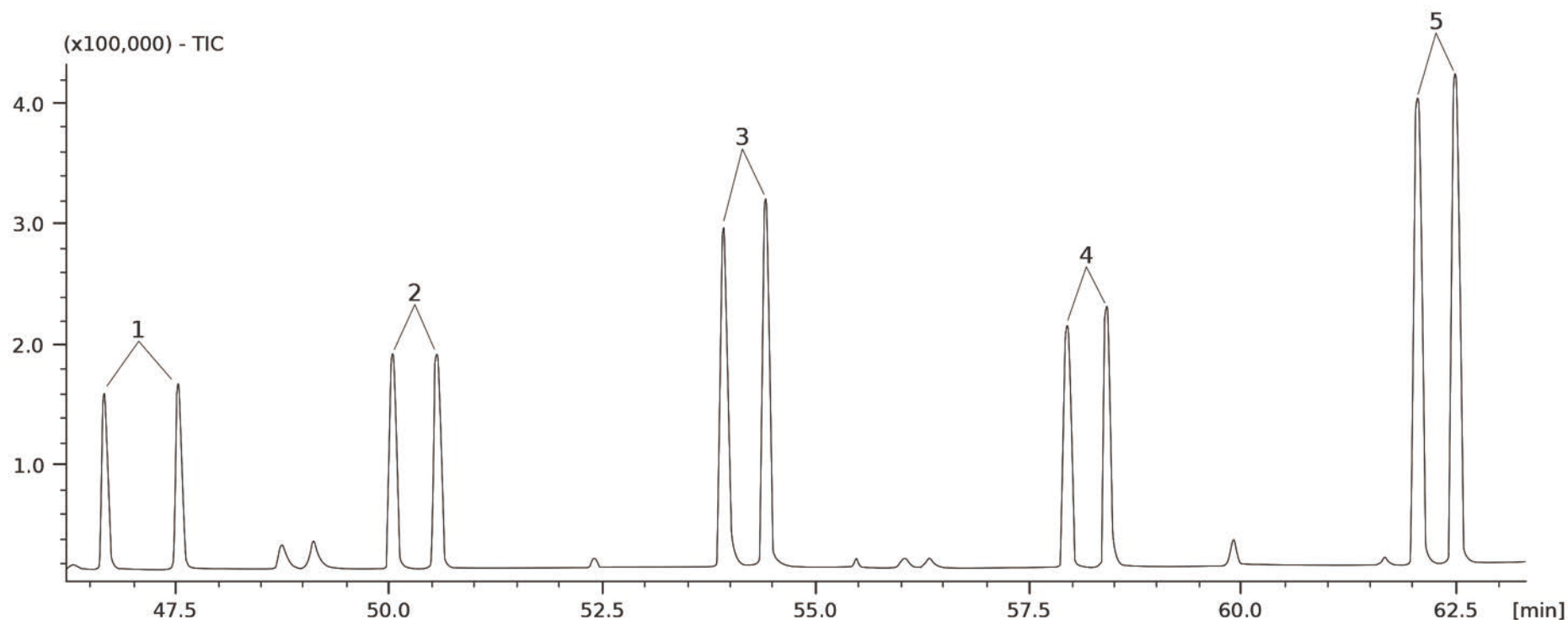

Acknowledgement: Prof. C.Bicchi, C.Cordero et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.