

异薄荷醇，新薄荷醇和 新异薄荷醇

Chiral Separation on MEGA-DEX Columns


条件:

进样口: 分流 250°C, 1:50 分流比, 0.5µL 进样量

制样: 标准混合, 1mL 正己烷中溶解 1µL 样品

Column: **MEGA-DEX DAC Beta - 0.25mm, 0.25µm, 25m**

Catalog Code: DEX-DAC-B-025-025-25

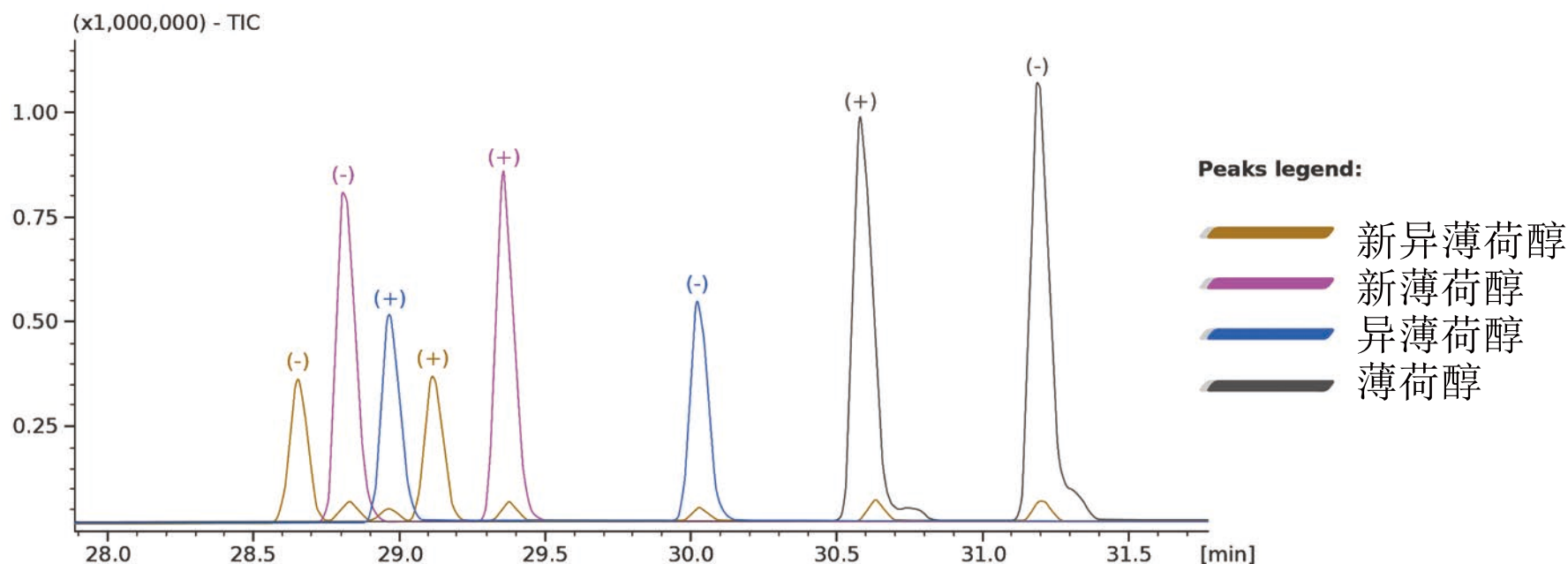

Acknowledgement: Prof. C.Bicchi, C.Cordero et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.