


Class 2 Mix B standard solution

1. Hexane
2. Nitromethane
3. Chloroform
4. 1,2-Dimethoxyethane
5. Trichloroethene
6. Pyridine
7. 2-Hexanone
8. Tetralin

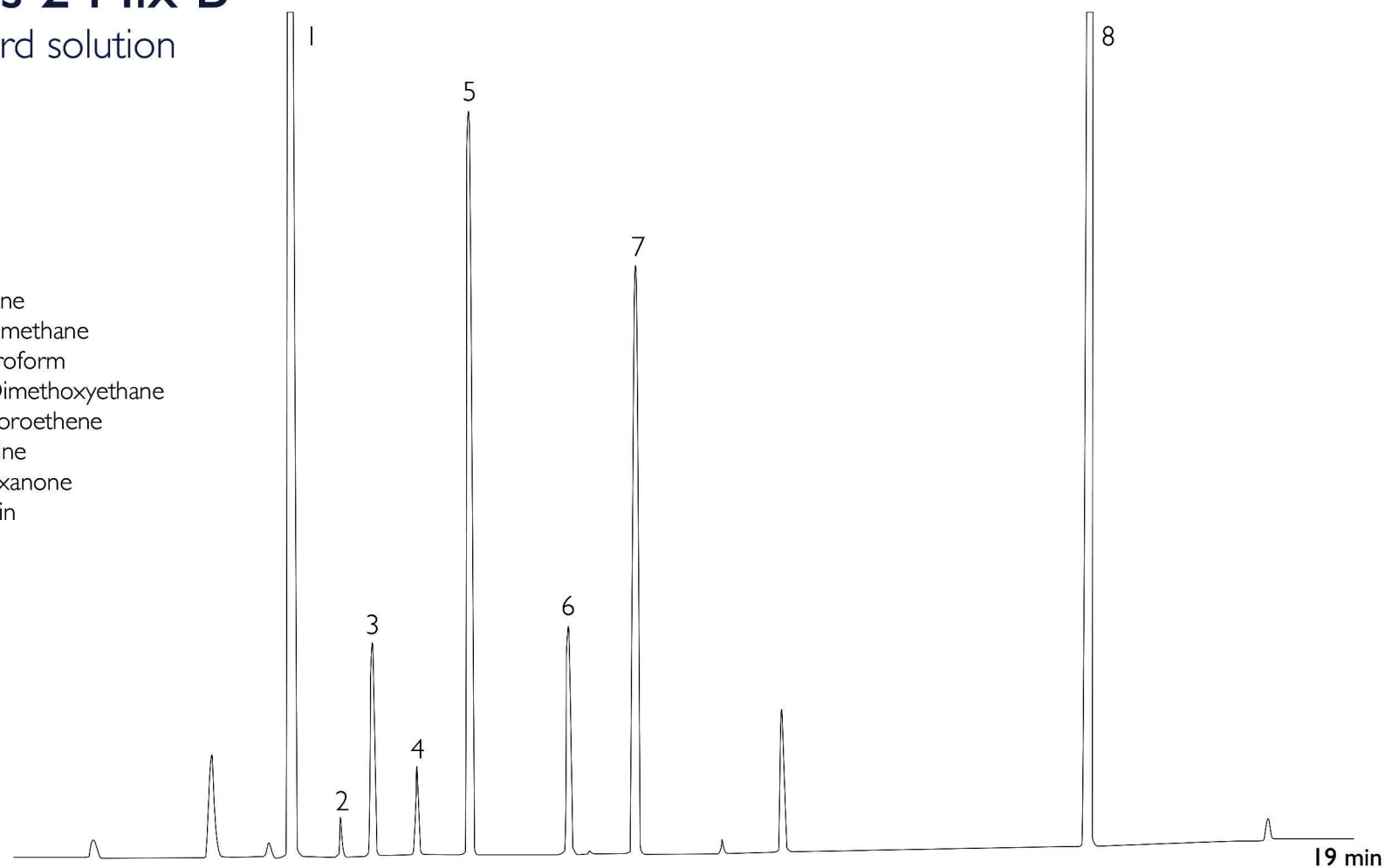

Analysis carried out with DANI Master-SHS and Master-GC system.

Column: **MEGA-624 - 0.18mm, 1.00 μ m, 20m**

Catalog Code: C-624-018-100-20

(USP G43 stationary phase)

条件:

进样口: 顶空平衡温度为80 $^{\circ}$ C, 1mL进样量,

平衡时间为40分钟, 传送温度为90 $^{\circ}$ C

检测器: FID 200 $^{\circ}$ C

柱温箱程序: 35 $^{\circ}$ C, 10 $^{\circ}$ C/min

载气: 氦气, 100kPa

Sample:

Residual Solvents - Class 2 Mix B, headspace sampling.