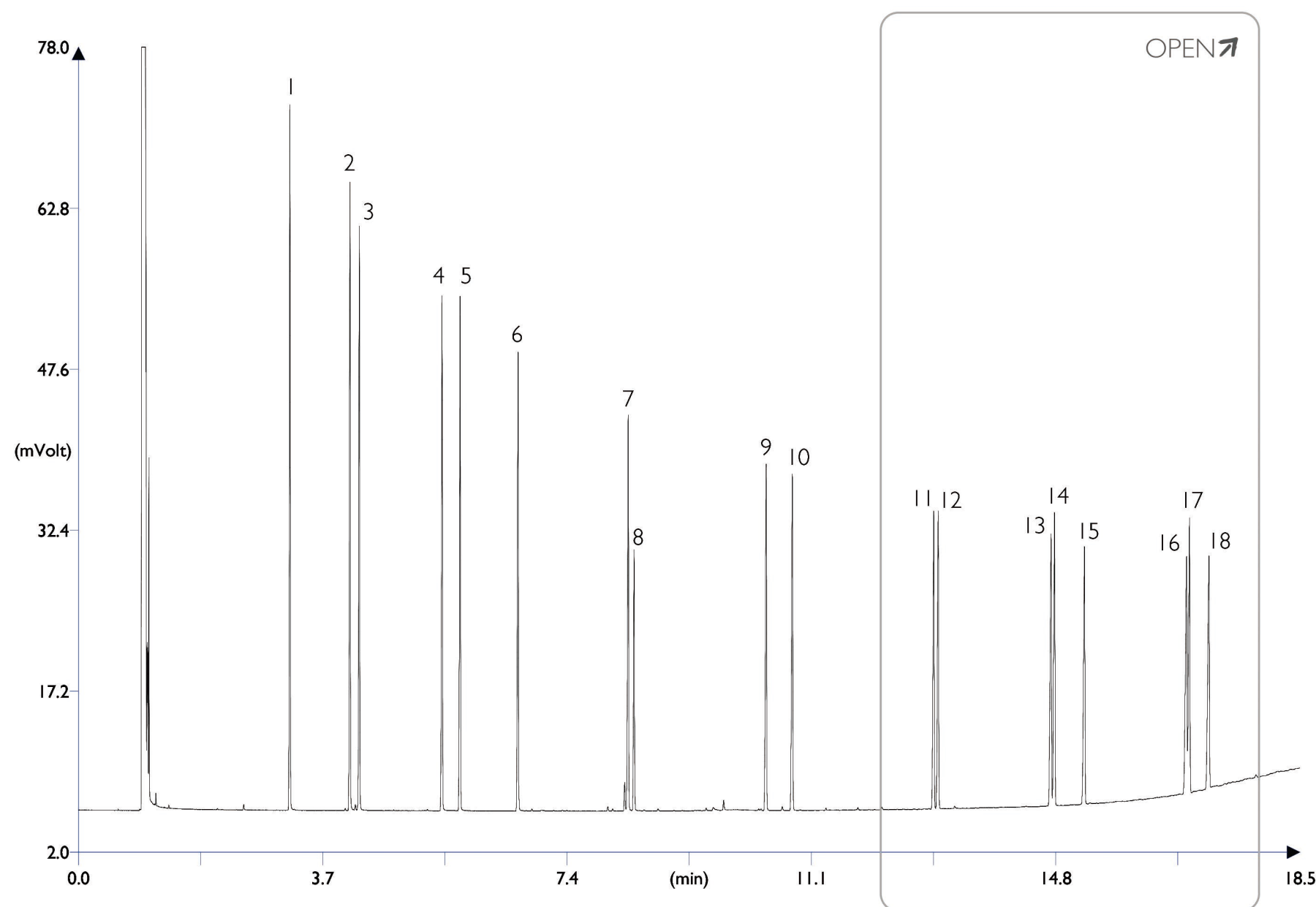




Column: **MEGA-5 MS Xil** - 0.15mm, 0.10µm, 30m

Catalog Code: MS-5XIL-015-010-30

条件:

进样口: 分流 230°C, 200mL/min 分流比, 0.2µL 进样量

检测器: FID 330°C

柱温箱程序: 90°C-300°C, 0.5°C/min

载气: 氢气, 240kPa

Sample:

8270 Calibration Mix #5, Revised 2000µg/mL in Methylene Chloride.

EPA method 610, 8100, 8270

- | | |
|------------------------|----------------------------|
| 1. Naphthalene | 10. Pyrene |
| 2. 1-methylnaphthalene | 11. Benzo[a]anthracene |
| 3. 2-methylnaphthalene | 12. Chrysene |
| 4. Acenaphthylene | 13. Benzo[b]fluoranthene |
| 5. Acenaphthene | 14. Benzo[k]fluoranthene |
| 6. Fluorene | 15. Benzo[a]pyrene |
| 7. Phenanthrene | 16. Indeno[1,2,3-cd]pyrene |
| 8. Anthracene | 17. Dibenzo[a,h]anthracene |
| 9. Fluoranthene | 18. Benzo[g,h,i]perylene |