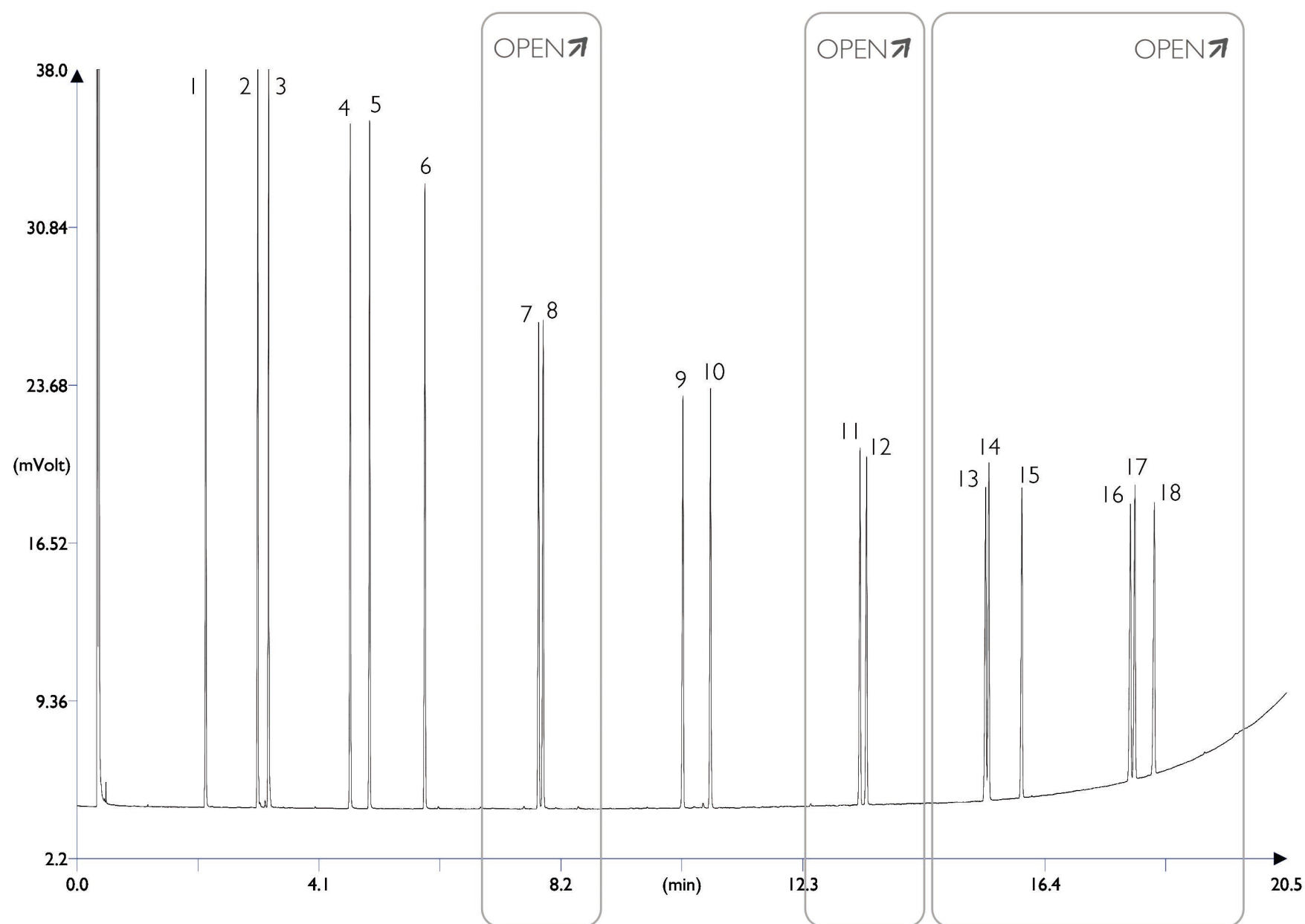




Column: **MEGA-PAH - 0.15mm, 0.15 μ m, 15m**

Catalog Code: C-PAH-015-015-15

条件:

进样口: 分流 230 $^{\circ}$ C, 300mL/min 分流比, 0.2 μ L 进样量

检测器: FID 330 $^{\circ}$ C

柱温箱程序: 90 $^{\circ}$ C-340 $^{\circ}$ C, 12 $^{\circ}$ C/min

载气: 氢气, 170kPa

Sample:

8270 Calibration Mix #5, Revised 2000 μ g/mL in Methylene Chloride.

EPA method 610, 8100, 8270

- | | |
|------------------------|----------------------------|
| 1. Naphthalene | 10. Pyrene |
| 2. 1-methylnaphthalene | 11. Benzo[a]anthracene |
| 3. 2-methylnaphthalene | 12. Chrysene |
| 4. Acenaphthylene | 13. Benzo[b]fluoranthene |
| 5. Acenaphthene | 14. Benzo[k]fluoranthene |
| 6. Fluorene | 15. Benzo[a]pyrene |
| 7. Phenanthrene | 16. Indeno[1,2,3-cd]pyrene |
| 8. Anthracene | 17. Dibenzo[a,h]anthracene |
| 9. Fluoranthene | 18. Benzo[g,h,i]perylene |