

Column: **MEGA-PAH FAST - 0.10mm, 0.10 μ m, 15m**

Catalog Code: F-PAH-010-010-15

EPA method 8100 (PAHs)


条件:

进样口: 分流 280 $^{\circ}$ C

检测器: FID 310 $^{\circ}$ C

柱温箱程序: 120 $^{\circ}$ C-30 $^{\circ}$ C/min-220 $^{\circ}$ C-15 $^{\circ}$ C/min-300 $^{\circ}$ C,

载气: 氢气, 0.6mL/min

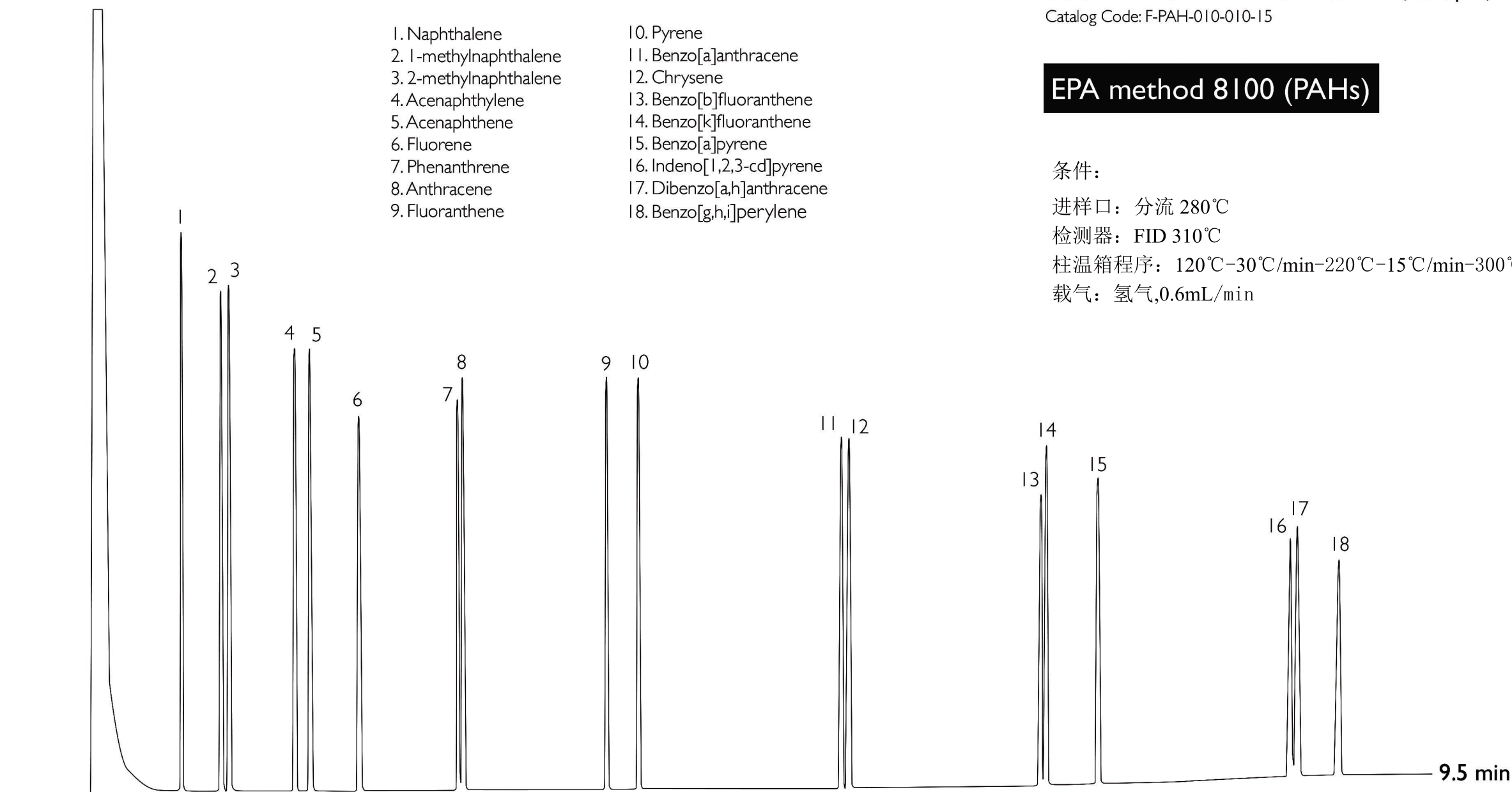

Analysis carried out with DANI Master GC.