


二甲苯异构体 (99%对二甲苯)

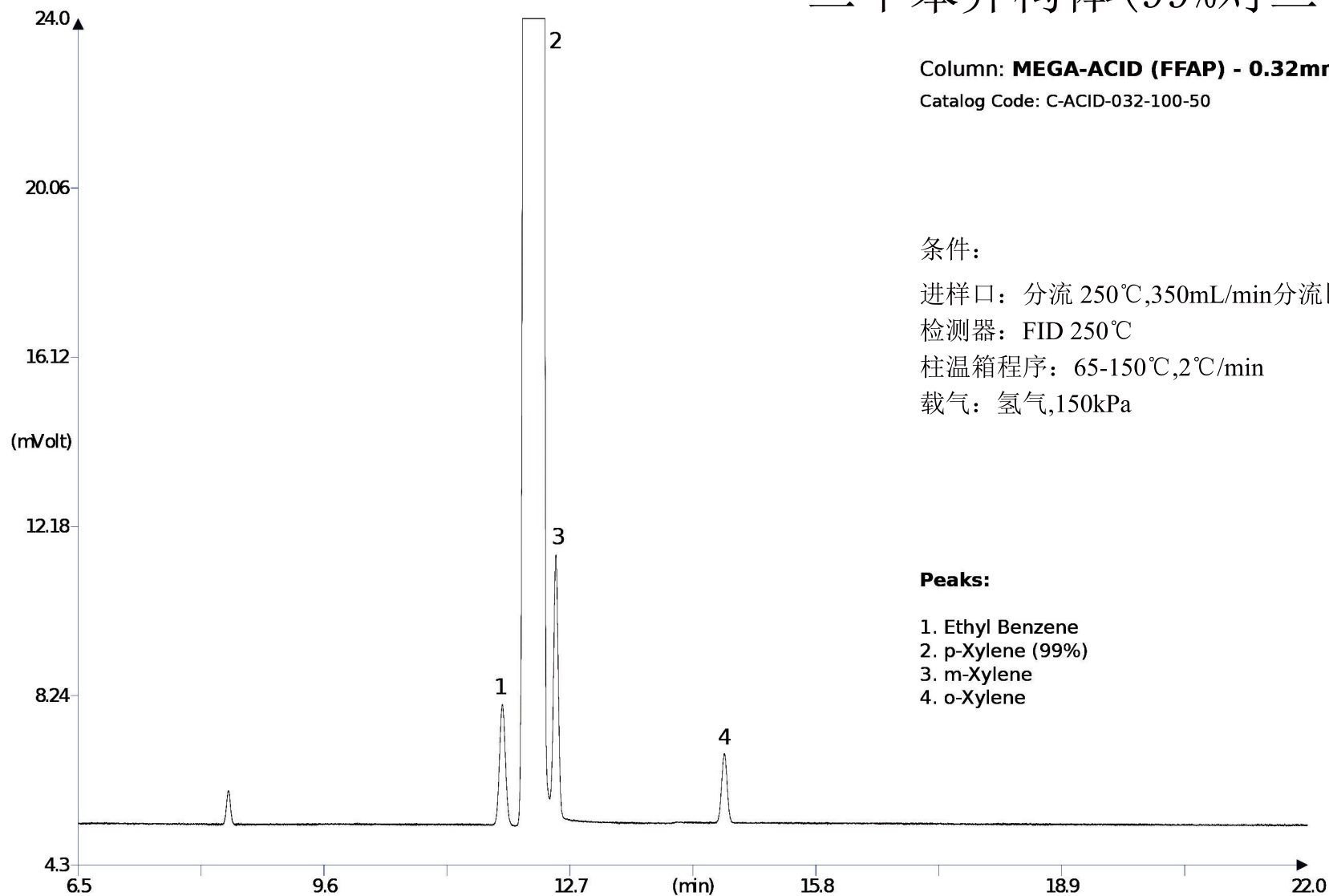

Column: **MEGA-ACID (FFAP) - 0.32mm, 1.00µm, 50m**

Catalog Code: C-ACID-032-100-50

条件:

进样口: 分流 250°C, 350mL/min 分流比, 0.5µL 进样量

检测器: FID 250°C

柱温箱程序: 65-150°C, 2°C/min

载气: 氢气, 150kPa

Peaks:

1. Ethyl Benzene
2. p-Xylene (99%)
3. m-Xylene
4. o-Xylene