

条件:

进样口: 分流 250℃, 1:50分流比, 0.5μL进样量

检测器: FID 250℃


制样: 标准混合, 1mg/mL溶解在正己烷

沉香醇, 芳樟氧化物

Chiral Separation on MEGA-DEX Columns

Column: **MEGA-DEX DMT Beta - 0.25mm, 0.25μm, 25m**

Catalog Code: DEX-DMT-B-025-025-25

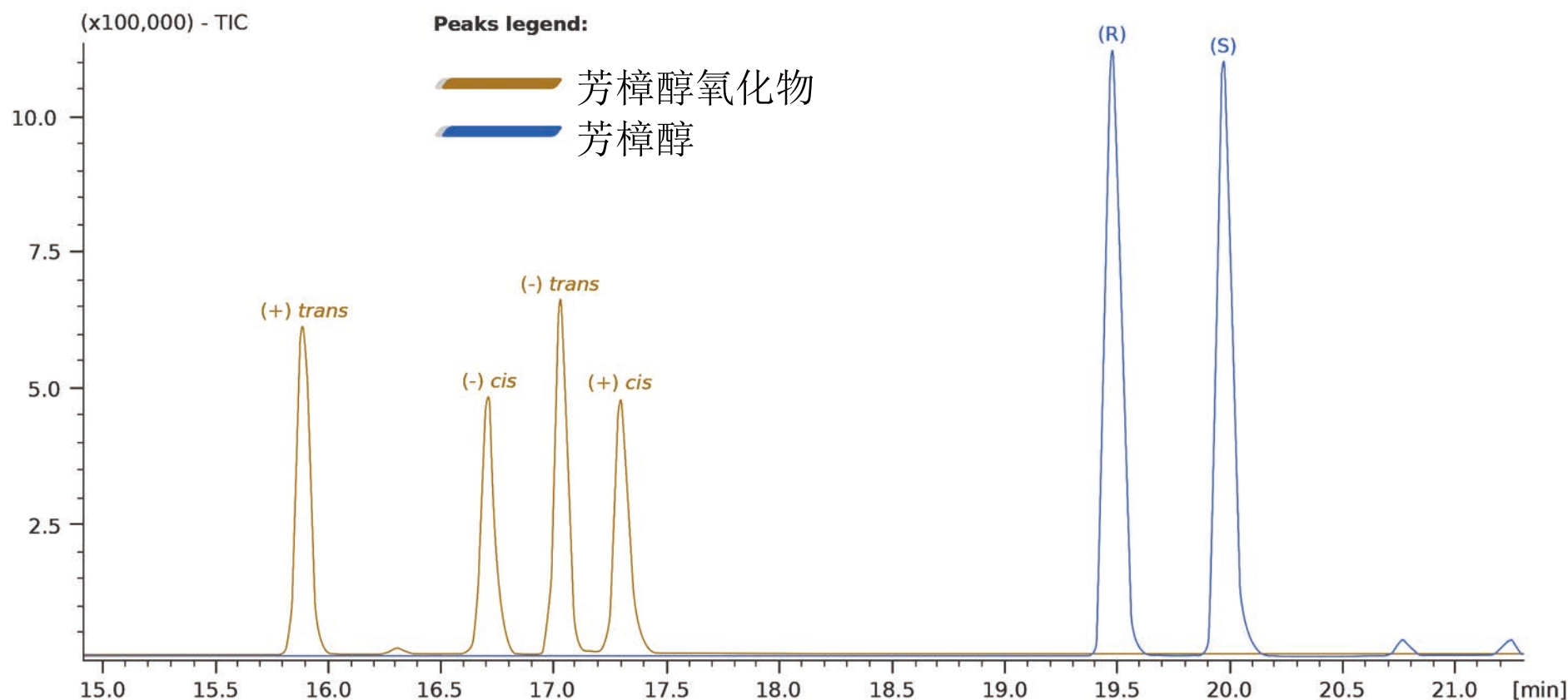

Acknowledgement: Prof. C.Bicchi, C.Cordero et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.