


咖啡香气中的硫化物

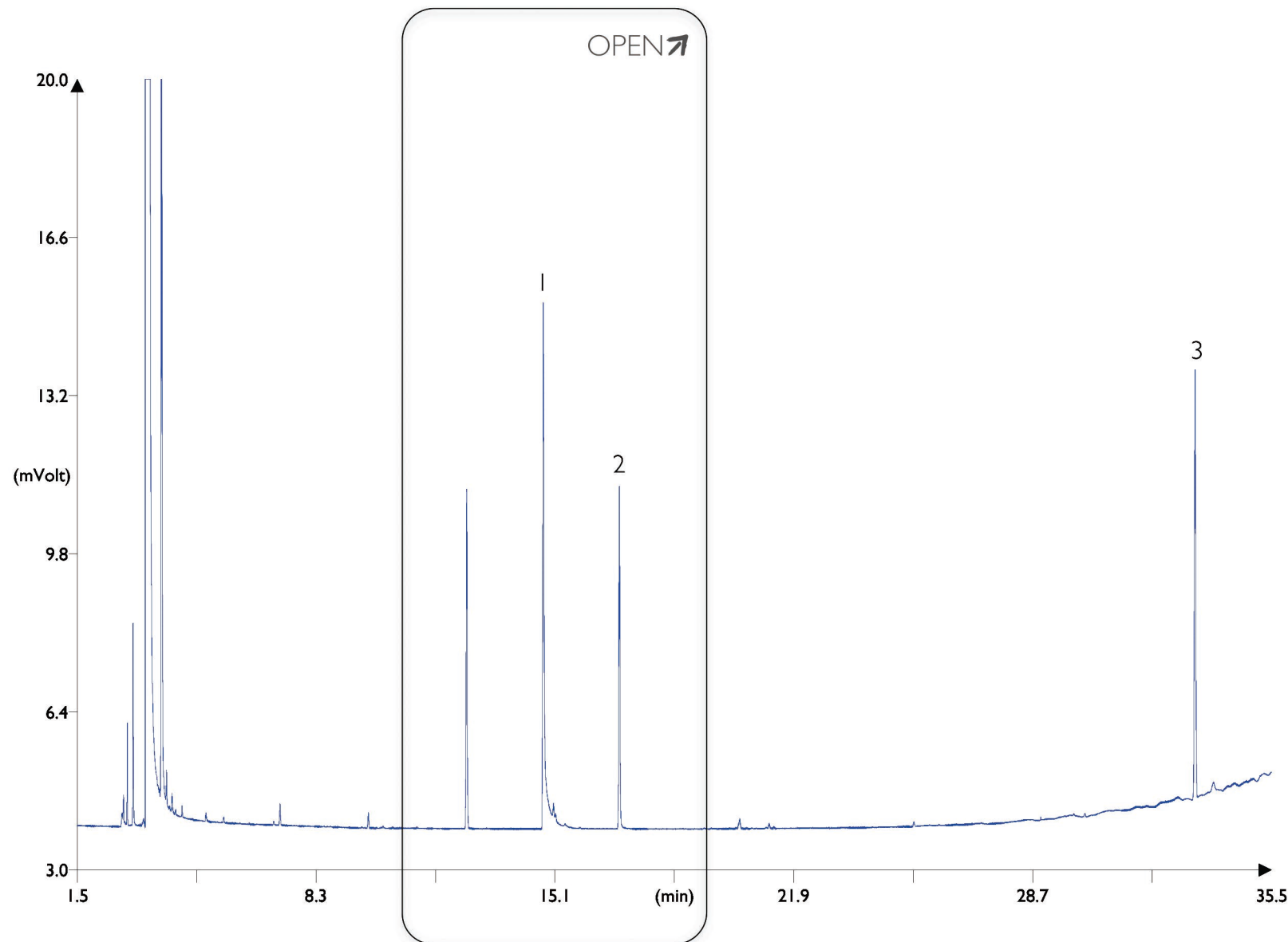

Column: **MEGA-ACID (FFAP)** - 0.25mm, 0.25 μ m, 60m

Catalog Code: C-ACID-025-025-60

USP G25, G35, G39 stationary phase

条件:

进样口: 分流 250 $^{\circ}$ C, 1:15分流比, 1.5 μ L进样量

检测器: FID 250 $^{\circ}$ C

柱温箱程序: 60-250 $^{\circ}$ C, 5 $^{\circ}$ C/min

载气: 氢气, 120kPa

Sample:

100ppm each compound in Acetone.

1. 2-Furfurylthiol
2. 3-Mercapto-3-methylbutyl formate
3. Furfuryl sulfide