

Column: **MEGA-LAP - 0.32mm, 0.10 μ m, 25m**

Catalog Code: C-LAP-032-010-25

条件:

进样口: 分流 300 $^{\circ}$ C

检测器: FID 300 $^{\circ}$ C


柱温箱程序: 220-280 $^{\circ}$ C, 5 $^{\circ}$ C/min

载气: 氢气, 2mL/min

Sample:

TMS derivatives of extracted from Arabica roasted coffee.

1. Caffeol
2. Cafestol
3. 16-O-methylcafestol
4. α -Cholesterol (I.S.)
5. 22,23-dihydro-brassicasterol
6. Campesterol
7. Stigmasterol
8. Sitosterol
9. Δ (5)-avenasterol
10. Δ (7)-stigmastenol
11. Δ (7)-avenasterol

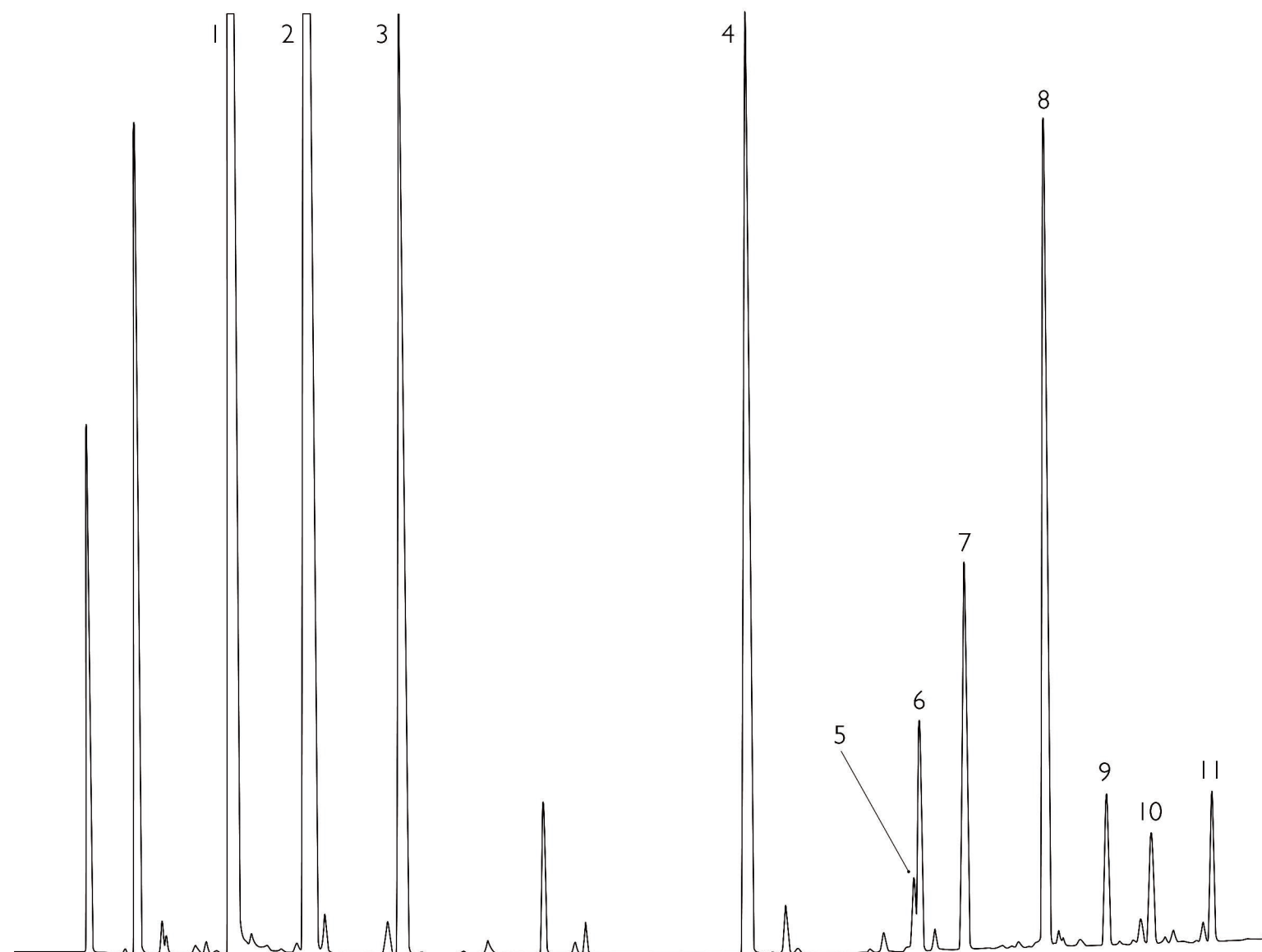

Acknowledgement: **Prof. C. Mariani**, Stazione Sperimentale Oli e Grassi, Via Giuseppe Colombo, 79 - 20133 - Milano - Italy.