

茨酮(樟脑)

Chiral Separation on MEGA-DEX Columns

Phase: **MEGA-DEX DET Beta**

Catalog Code: DEX-DET-B-025-025-25

Phase: **MEGA-DEX DMT Beta**

Catalog Code: DEX-DMT-B-025-025-25

Phase: **MEGA-DEX DAC Beta**

Catalog Code: DEX-DAC-B-025-025-25


Columns Dimensions: **0.25mm, 0.25 μ m, 25m.**

条件:

进样口: 分流250 $^{\circ}$ C, 0.5 μ L, 1:50分流比

检测器: FID 250 $^{\circ}$ C

制样: 标准混合, 1mg的样品溶解在1mL的正己烷

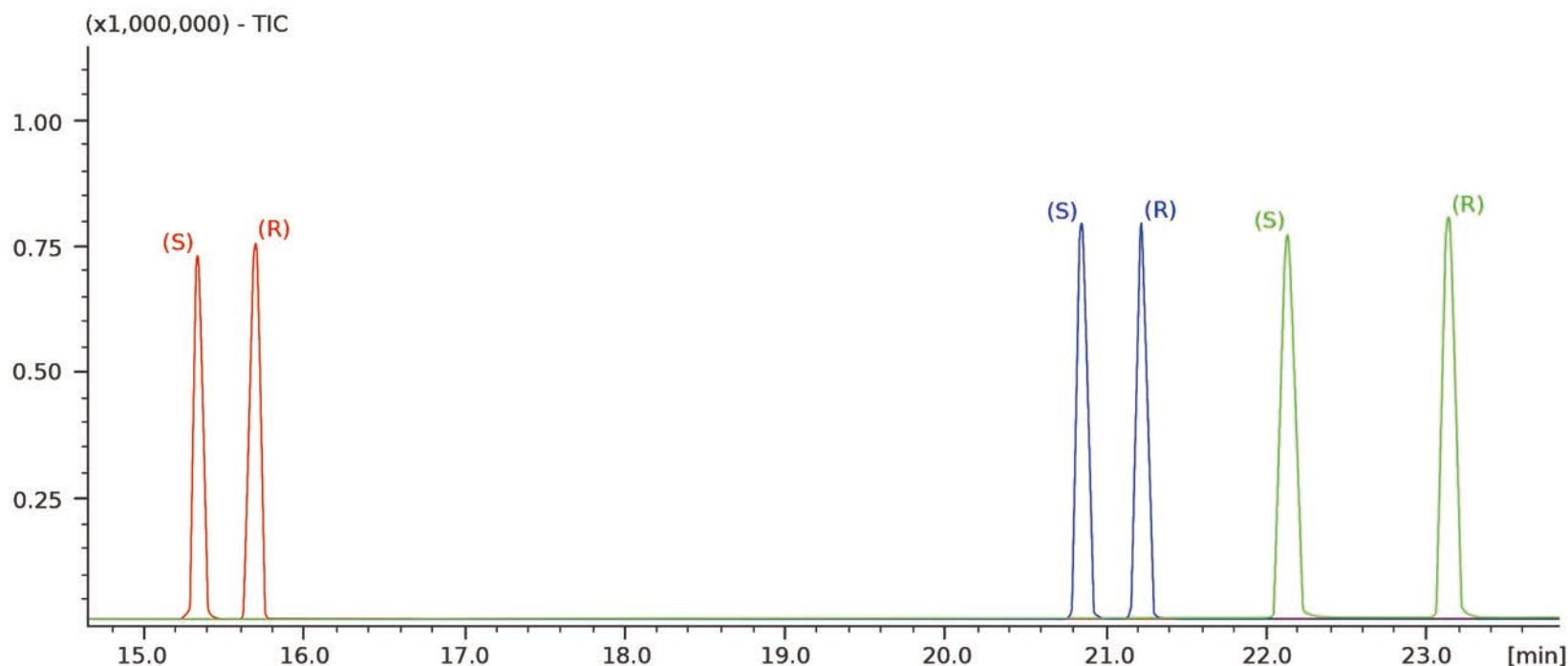

Acknowledgement: Prof. C.Bicchi, C.Cordero et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.