

卤代水化合物

Column: **MEGA-DAI 1 - 0.32mm, 5.00µm, 30m**

Catalog Code: C-DAI1-032-500-30

MEGA-DAI Columns (Direct Aqueous Injections) allow direct introduction of aqueous solutions into the column without sample preparation.


条件:

进样口: 冷柱进样

检测器: ECD 300°C, 尾吹为50mL/min, 氮气

柱箱程序: 恒温104°C

载气: 氦气, 0.8mL/min

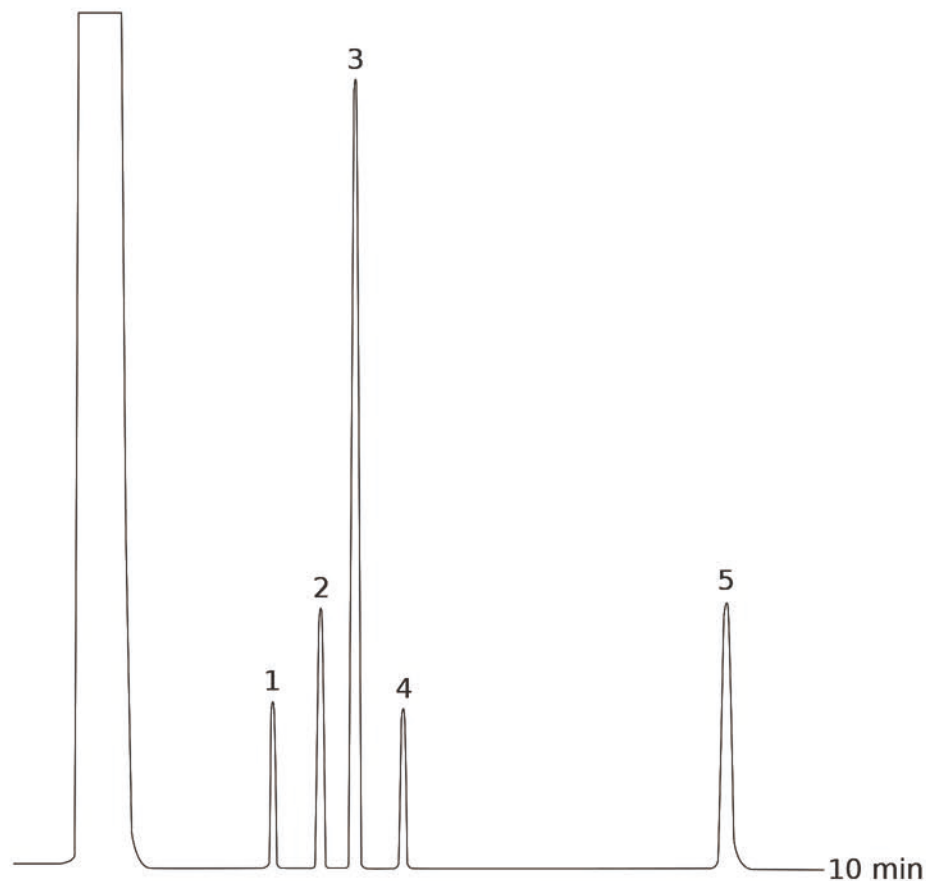

峰:

1. 氯仿
2. 甲基氯仿
3. 四氯化碳
4. 三氯乙烯
5. 四氯乙烯