


氯苄苯丙胺手性分离

dex xeb
chiral columns

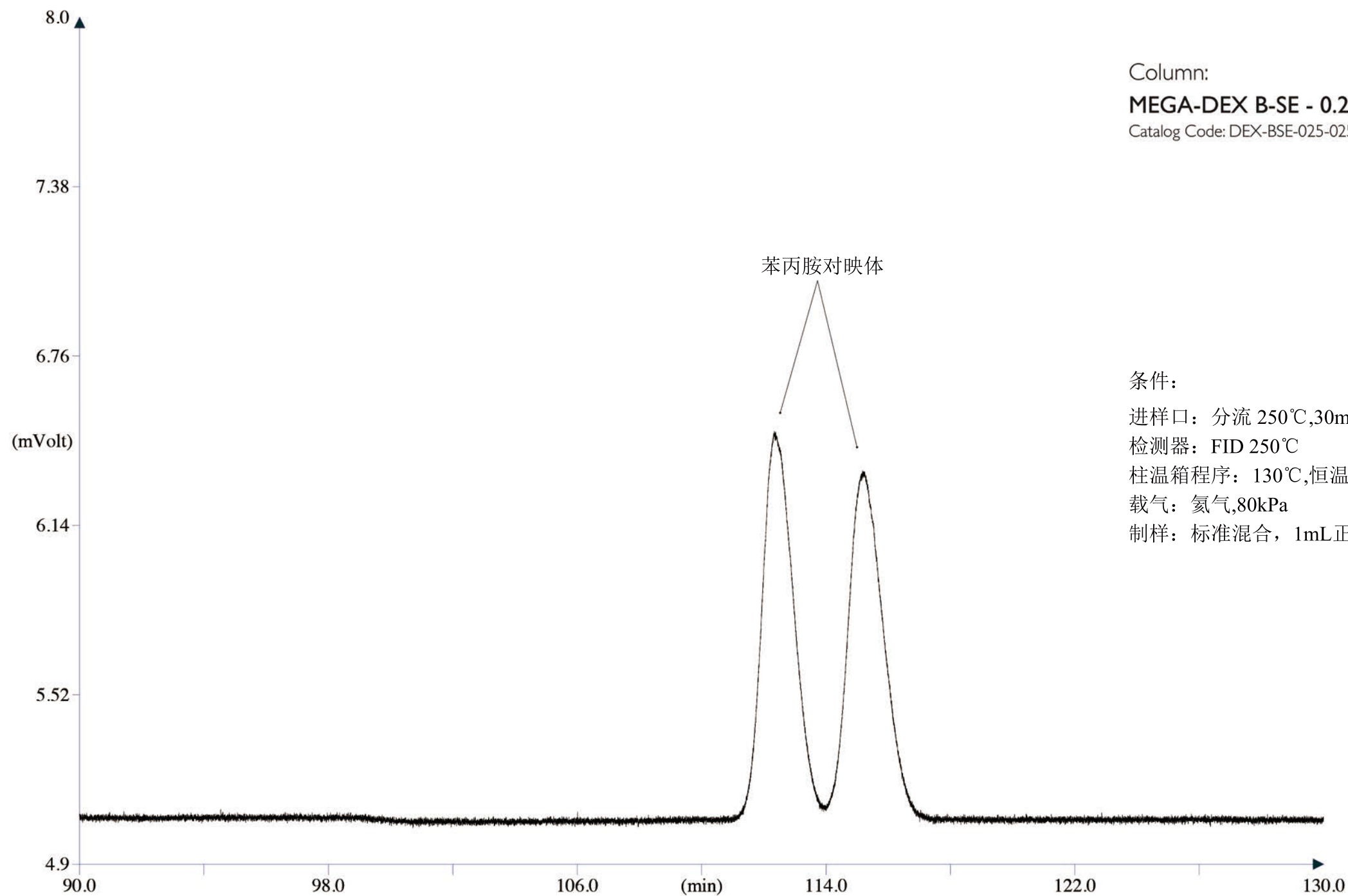

Column:

MEGA-DEX B-SE - 0.25mm, 0.25 μ m, 15m

Catalog Code: DEX-BSE-025-025-15

条件:

进样口: 分流 250°C, 30mL/min 分流比, 1 μ L 进样量

检测器: FID 250°C

柱温箱程序: 130°C, 恒温模式

载气: 氮气, 80kPa

制样: 标准混合, 1mL 正己烷中溶解 1 μ L 样品

MEGA®

www.mega.mi.it