

酞酸酯

EPA Method 606

Column: **MEGA-1 - 0.53mm, 1.50 μ m, 15m**

Catalog Code: S-1-053-150-15

条件:


进样口: 分流 250 $^{\circ}$ C, 1.0 μ L进样量

检测器: FID 300 $^{\circ}$ C

柱温箱程序: 150-270 $^{\circ}$ C, 15 $^{\circ}$ C/min

载气: 氦气, 20mL/min

制样: 2000 μ g/mL在甲醇

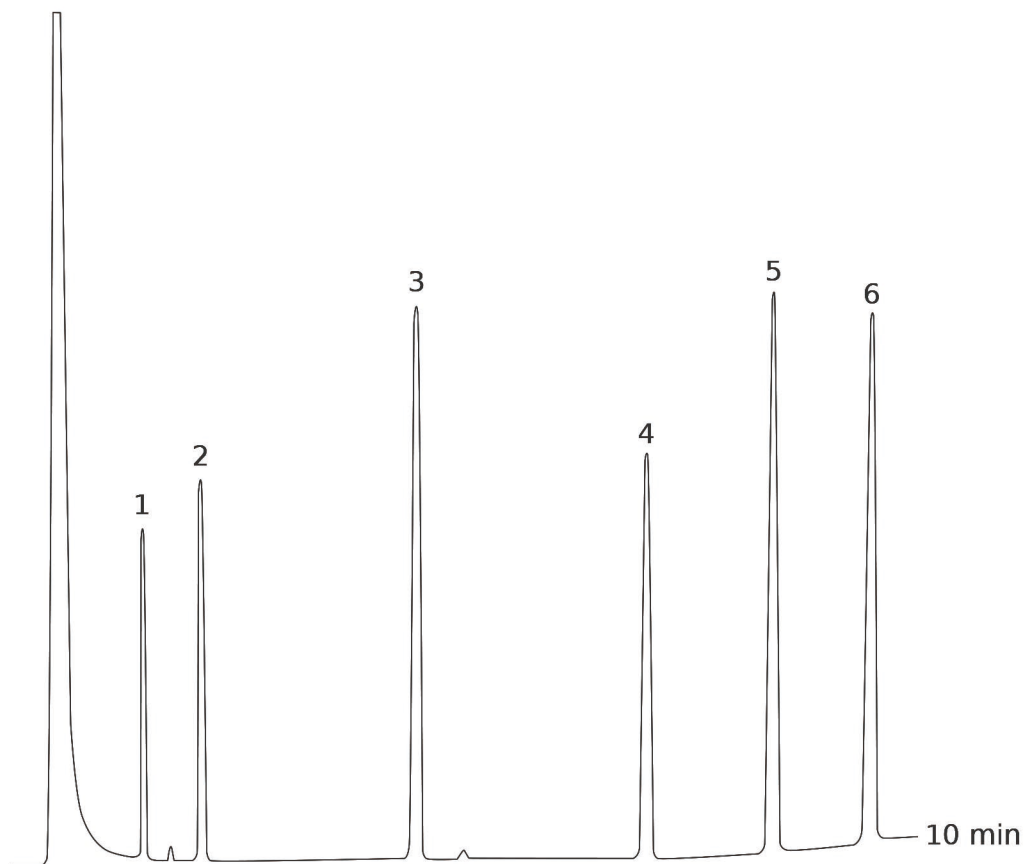

Peaks:

1. Dimethyl phthalate
2. Diethyl phthalate
3. Di-n-butyl phthalate
4. Butyl benzyl phthalate
5. Bis (2-ethylhexyl) phthalate
6. Di-n-octyl phthalate