


Column: **MEGA-624 - 0.18mm, 1.00µm, 20m**
Catalog Code: C-624-018-100-20

条件:
进样口: 不分流 250℃,吹扫补集
检测器: TOF MS
柱温箱程序: 35-240℃,20℃/min
载气: 氦气, 0.8mL/min

EPA method 524.3

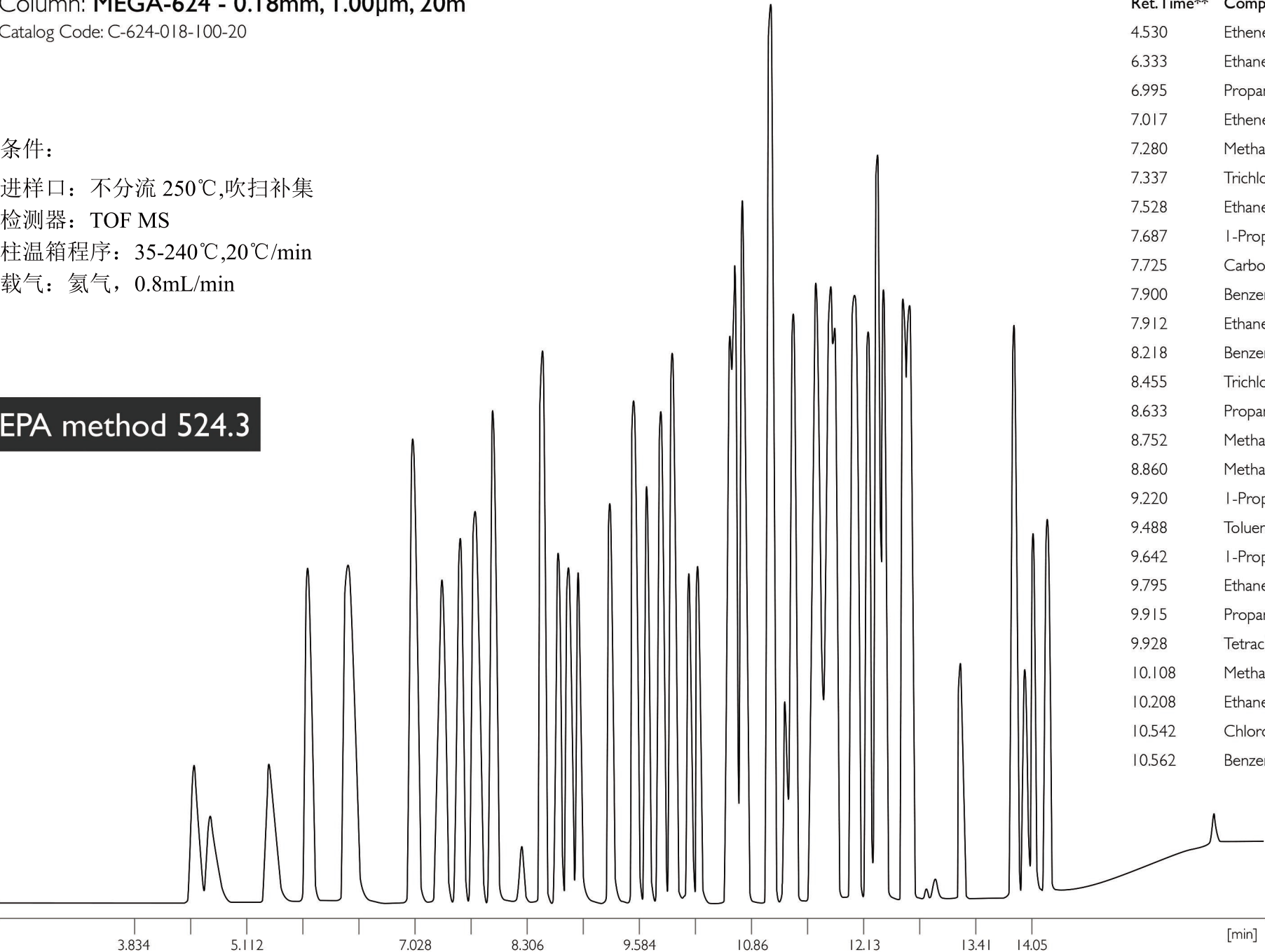

Ret.Time**	Compound
4.530	Ethene, 1,1-dichloro-
6.333	Ethane, 1,1-dichloro-
6.995	Propane, 2,2-dichloro-
7.017	Ethene, 1,2-dichloro-
7.280	Methane, bromochloro-
7.337	Trichloromethane
7.528	Ethane, 1,1,1-trichloro-
7.687	1-Propene, 1,1-dichloro-
7.725	Carbon Tetrachloride
7.900	Benzene
7.912	Ethane, 1,2-dichloro-
8.218	Benzene, 1,4-difluoro- (I.S. 1)
8.455	Trichloroethylene
8.633	Propane, 1,2-dichloro-
8.752	Methane, dibromo-
8.860	Methane, bromodichloro-
9.220	1-Propene, 1,3-dichloro-, (Z)-
9.488	Toluene
9.642	1-Propene, 1,3-dichloro-, (E)-
9.795	Ethane, 1,1,2-trichloro-
9.915	Propane, 1,3-dichloro-
9.928	Tetrachloroethylene
10.108	Methane, dibromochloro-
10.208	Ethane, 1,2-dibromo-
10.542	Chlorobenzene-d5 (I.S. 2)
10.562	Benzene, chloro-

Ret.Time**	Compound
10.628	Ethylbenzene
10.713	p-Xylene
11.013	o-Xylene
11.027	Styrene
11.278	Benzene, 1-methylethyl-
11.503	Ethane, 1,1,2,2-tetrachloro-
11.538	Benzene, bromo-
11.547	Propane, 1,2,3-trichloro-
11.583	Benzene, propyl-
11.668	Benzene, 1-chloro-2-methyl-
11.707	Benzene, 1,2,4-trimethyl-
11.742	Benzene, 1-chloro-3-methyl-
11.955	Benzene, tert-butyl-
11.990	Benzene, 1,3,5-trimethyl-
12.117	Benzene, 1-methylpropyl-
12.217	Benzene, 1-methyl-4-(1-methylethyl)-
12.227	Benzene, 1,4-dichloro-
12.273	1,4-Dichlorobenzene-d4 (I.S. 3)
12.290	Benzene, 1,3-dichloro-
12.525	Benzene, butyl-
12.577	Benzene, 1,2-dichloro-
13.165	Propane, 1,2-dibromo-3-chloro-
13.768	Benzene, 1,2,3-trichloro-
13.885	1,3-Butadiene, 1,1,2,3,4,4-hexachloro-
13.975	Naphthalene
14.140	Benzene, 1,2,4-trichloro-

**:[min].

*: analysis carried out with DANI Master GC, DHS and TOF-MS system.
Acknowledgment: Alessandra Mantegazza, Daniela Cavagnino,
courtesy of DANI Instruments S.p.A., Cologno Monzese (MI), Italy.